

# Drugs of Abuse Rapid Drug Tests

**Neopharma**  
Healthcare & Innovation

- ▶ Quick results
- ▶ Simple and easy to use
- ▶ Sensitive
- ▶ No cross-reactivity
- ▶ Available for up to 14 different drugs of abuse under Neotest brand



Neotest's Rapid Drug Tests for Drugs of Abuse are quick, sensitive, accurate, and reliable. These rapid tests use immuno-chromatographic assay to provide information about the possible use of drugs of abuse. The product is manufactured by **Artron Laboratories Inc, Canada**.

Our Company offers Rapid Drug Tests for the following Drugs of Abuse:

Amphetamines (AMP), Barbiturates (BAR), Benzodiazepines (BZO), Benzoylcegonine (BCOC), Cocaine (COC), Ecstasy (MDMA), Methamphetamines (MET), Opiates (OPI), Buprenorphine (BUP), Oxycodone (OXY), Tricyclic antidepressants (TCA), Marijuana (THC), Morphine (MOR) Ketamine (KET) and Methadone (MTD)


## Available in 4 Formats




**Stick**



**Dip Card**

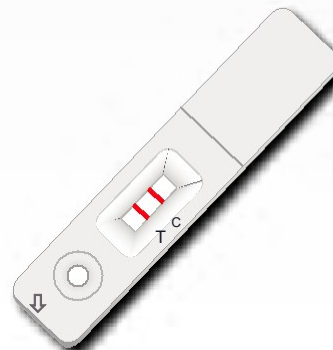

**Cassette**

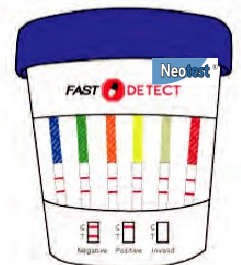

**Cup**

NEOPHARMA BIOTECH ASIA SDN BHD. (648086-W)

Unit 3.1, Level 3, Block A, Dataran PHB Persiaran Golf, Saujana Resort, Seksyen U2, 40150 Shah Alam, Selangor Darul Ehsan.

Tel: +603-7842 6787 (Hunting) Fax: +603-7842 5787 [www.neopharmaasia.com](http://www.neopharmaasia.com)

## How it Works

Our Rapid Drug Tests are manufactured in conformity with “cut-off” levels established by The Substance Abuse and Mental Health Services Agency (SAMHSA). If the level of drugs in the sample is above the cut-off level, the test will yield a positive result, and as long as the level of drug in the sample is below the cut-off level, the test will yield a negative result.

### Drug Clearance Rates

Drug	Time until first detection	Max. time drug is detectable
Amphetamine	4 - 6 hours	2 - 3 days
Barbituates	2 - 4 hours	2 - 10 days
Bezoylecgonine	2 - 5 hours	3 - 4 days
Buprenorphine	1 - 4 hours	1 - 3 days
Benzodiazepine	2 - 7 hours	1 - 4 days
Cocaine	2 - 6 hours	2 - 3 days
Ecstasy	2 - 7 hours	2 - 4 days
Methamphetamine	4 - 6 hours	2 - 3 days
Morphine	1 - 2 hours	8 - 24 hours
Methadone	3 - 8 hours	1 - 3 days
Opiates	2 - 6 hours	1 - 3 days
Oxycodone	4 - 8 hours	3 - 4 days
Ketamine	4 - 8 hours	2 - 8 days
Tricyclic Antidepressants	8 - 12 hours	2 - 7 days
Marijuana	1 - 3 hours	1 - 7 days

### Cut-Off Levels

Identifier	Drug Name	Cut-Off Level
AMP	Amphetamine	1000 ng/mL
BAR	Barbituates	300 ng/mL
BCOC	Bezoylecgonine	300 ng/mL
BUP	Buprenorphine	50 ng/mL
BZO	Benzodiazepine	300 ng/mL
COC	Cocaine	300 ng/mL
MDMA	Ecstasy	500 ng/mL
MET	Methamphetamine	1000 ng/mL
MOR	Morphine	300 ng/mL
MTD	Methadone	300 ng/mL
OPI	Opiates	2000 ng/mL
OXY	Oxycodone	1000 ng/mL
KET	Ketamine	1000ng/mL
TCA	Tricyclic Antidepressants	1000 ng/mL
THC	Marijuana	50 ng/mL

The drug clearance rate indicates how soon an individual may receive a positive test result after taking a particular drug. It also indicates how long the individual may continue to test positive after the last time the individual took the drug. However, these are only guidelines and times may vary significantly from these estimates based on how long the individual has been taking the drug, the amount of drugs used, and the individual’s metabolism.

## Quick Results

Neotest’s Rapid Drug Tests give quick results!

Results are easy to interpret.

One coloured band indicates a positive result, and two-coloured bands indicate a negative result.



Positive



Negative